

Polymer Chemie

Viscometry

Universität Stuttgart

Author 1:	Laurens Viehoff
Author 2:	Felix Roth
E-Mail 1:	st183688@stud.uni-stuttgart.de
E-Mail 2:	st188157@stud.uni-stuttgart.de
Student-ID 1:	3676761
Student-ID 2:	3711192
Assistentin:	Yanzhe Huang

Inhaltsverzeichnis

1	Theory	2
1.1	Viscosity	2
1.2	From viscosity to M_η	2
2	Experimental	4
3	Results	4
4	Questions	6
4.1	Question 1	6
4.2	Question 2	6
4.3	Question 3	6
4.4	Question 4	6

1 Theory

1.1 Viscosity

In a flowing liquid the brutto-fluid-flowing-peed is reversed proportional to the viscosity. This is because at the wall of the pipe or flask there is a flowgradient. This gradient can be described by many different fast flowing layers also shown in Scheme 1.

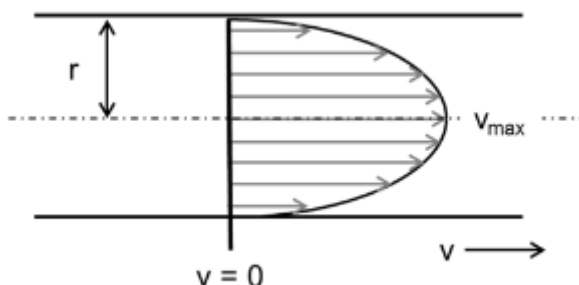


Abb. 1: The picture shows a liquid flöwing through a pipe with the simplifycation of the fluid layers. r - is the capillaryradius and v is the flowspeed [polyskript]

Viscosity in polymers is influenced by 5 key factors. The first being the chain lenght of the Polymer, the longer the chains the more viscouse the fluid gets. This is due to the Intramolecular forces.

The second influence is the chemical structure of the polymer, different functional groups mean different intermolecular forces. For example if the side chain has a substituent which can make H-Bonds the viscosity gets higher. In total it can be stated, the stronger the intramolecular forces the higher the viscosity.

The third one is the influence of the solvent. Because different solvents have different properties. If the solvent is better the tangle gets more swallen and the viscosity gets higher.

The fourth is the influence of the concentration. This means with different concentration the viscosity changes. This is due to higher concentration leads to more intermolecular forces which raises the viscosity.

The last but not least influence is the temperature. Because the viscosity changes with different temperatures because the movement of the molecules changes. The temperature also destroys the intermolecular forces which leads to a less tangeled structure wich also leads to a less rigid polymer. This results in a lower viscosity.

1.2 From viscosity to M_η

The viscosity is defined after Hagen-Poiseuille by following equation:

$$\eta = \frac{\pi \cdot r^4 \cdot t \cdot \Delta p}{8 \cdot l \cdot V} \quad (1)$$

In which the Δp is defined as:

$$\Delta p = \rho \cdot g \cdot h \quad (2)$$

In these equations r is the capillary radius, t is the time the polymer needs to flow through the capillary, Δp is the pressure difference in the capillary, l is the length of the capillary, V is the Volume, ρ is the density of the solution and h is the height of the measured part.

After analyzing different polymers four different types of Viscosities could be determined. With c for concentration of the polymer solution.

$$\eta_{\text{rel}} = \frac{\eta}{\eta_0} \quad (3)$$

$$\eta_{\text{sp}} = \eta_{\text{rel}} - 1 \quad (4)$$

$$\eta_{\text{red}} = \frac{\eta_{\text{sp}}}{c} \quad (5)$$

$$\eta = \frac{\ln \eta_{\text{rel}}}{c} \quad (6)$$

If equation 1 is taken into account η_{sp} is changed to:

$$\eta_{\text{sp}} = \frac{t\rho - t_0\rho_0}{t_0\rho_0} \quad (7)$$

And because of the low concentration the density of the solution is almost equal to the density of the solvent. Therefore the equation can be simplified.

$$\eta_{\text{sp}} \approx \frac{t - t_0}{t_0} \quad (8)$$

To get the viscoussmiddeled molecular Volume the Staudinger-index is used, and it is connected via the Marc-Howink-Equation.

$$[\eta] = K \cdot \bar{M}_\eta^a \quad (9)$$

This equation is used in combination with the equatoin for the viscositynumber.

$$[\eta] = \lim_{c \rightarrow 0} (\eta_{\text{red}}) = \lim_{c \rightarrow 0} \left(\frac{\eta_{\text{sp}}}{c} \right) \quad (10)$$

So by extrapolating the results η_{red} can be determined and via this $[\eta]$ and from that $\bar{M}_\eta$ can be determined. This is also used in our case in the experiment.

2 Experimental

The five different polymer solutions were already prepared and ready to use. One pure Toluene solution and four polymer solutions with different concentrations. All of these five solutions were analyzed with a Viscosimeter.

3 Results

In order to determine the viscosity-weighted molecular weight M_η , the specific viscosity η_{sp} is first calculated. According to Equation (8), this depends on the transit time t . The transit times and the respective mean values for all samples are listed in Table 1.

Tab. 1: Stoffmengenanteile und Massenanteile der Copolymerisation

Concentration [$10^{-3} \frac{\text{g}}{\text{ml}}$]	t_1 [s]	t_2 [s]	t_3 [s]	t_4 [s]	t_5 [s]	Average transit time [s]
0	23.40	23.36	23.37	23.37	23.37	23.374
1.25	24.99	24.99	24.99	24.98	24.98	24.986
2.5	26.86	26.81	26.79	26.79	26.78	26.806
5	32.06	32.01	32.03	32.10	32.06	32.052
10	40.42	40.37	40.35	40.35	40.36	40,370

The transit time for the sample with a concentration of 0 is used as the standard time t_0 . This results in a specific viscosity of

$$\eta_{\text{sp}} = \frac{24.986 \text{ s} - 23.374 \text{ s}}{23.374 \text{ s}} = 68.966 \cdot 10^{-3}$$

for the sample with a concentration of $1.25 \cdot 10^{-3}$. The specific viscosity is then divided by the concentration c , yielding the reduced viscosity η_{red} . For the first sample, this results in a reduced viscosity of

$$\eta_{\text{red}} = \frac{68.966 \cdot 10^{-3}}{1.25 \cdot 10^{-3} \frac{\text{g}}{\text{ml}}} = 55.172 \frac{\text{ml}}{\text{g}}.$$

The specific and reduced viscosities for the other samples are listed in Table 2.

Tab. 2: Specific and reduced viscosity at different polymer concentrations.

Sampleconcentration [$10^{-3} \frac{\text{g}}{\text{ml}}$]	Specific viscosity [10^{-3}]	Reduced viscosity [$\frac{\text{ml}}{\text{g}}$]
1.25	68.966	55.172
2.5	146.830	58.732
5	371.267	74.253
10	727.133	72.713

Next, the reduced viscosity is plotted against concentration, a linear regression is performed and the results are extrapolated to $c = 0 \frac{\text{g}}{\text{ml}}$. The graph is shown in Figure 2.

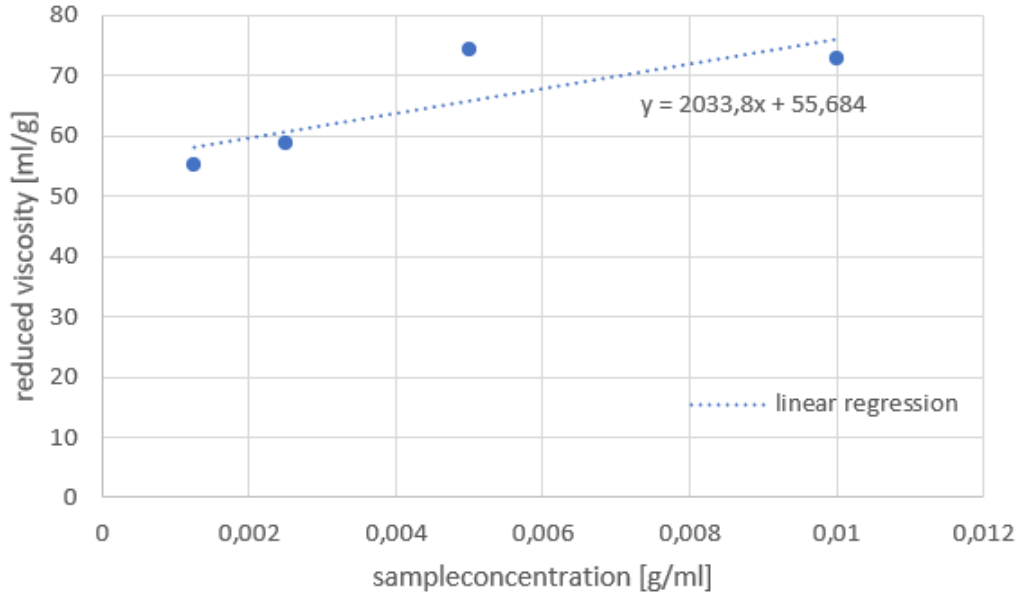


Abb. 2: Reduced viscosity η_{red} plotted against concentration c using a regression line $\eta_{\text{red}} = 2033.8 \frac{\text{ml}^2}{\text{g}^2} \cdot c + 55.684 \frac{\text{ml}}{\text{g}}$.

If a concentration of 0 is substituted into the linear equation, only the y-intercept remains, which in itself yields a viscosity number of $[\eta] = 55.684 \frac{\text{ml}}{\text{g}}$.

This quantity is proportional to the viscosity-weighted molecular weight according to Equation (9), which, when rearranged in terms of the viscosity-weighted molecular weight, yields

$$M_{\eta} = \sqrt[a]{\frac{[\eta]}{K}}. \quad (11)$$

Using the literature values of $K = 0.012 \frac{\text{ml}}{\text{g}}$ and $a = 0.71$, this results in a viscosity-weighted molecular weight of

$$M_{\eta} = \sqrt[0.71]{\frac{55.684 \frac{\text{ml}}{\text{g}}}{0.012 \frac{\text{ml}}{\text{g}}}} = 145933.088.$$

Therefore the viscosity-weighted molecular weight is $145\,933.088 \frac{\text{g}}{\text{mol}}$.

4 Questions

4.1 Question 1

In total it can be stated, all three of the given polymers have different properties. Polystyrene is a neutral polymer and behaves like a statistical tangle in a good solvent, the Staudinger-index only is a reference value. The Polystyrene-lithium is a living anionic polymer with Lithium as a counter ion. This Lithium-ion forms a polar Li-C-Bond in polar solvents, these are also called Units. These Units have a direct influence on the Staudinger-Index, so the Index gets higher. In the last case of the sulfonated Polystyrene is a bit more complex. PSS-Na is a Polyelectrolyte, therefore the SO_3^- is distributed over the chain. A Uniqueness is if this Polymer is put in demineralized water the negative charges repel each other. Therefore the polymer swells and the tangle gets bigger, the Staudinger index gets very high and is the highest from these three examples.

4.2 Question 2

At a certain concentration point things like the intramolecular forces, the hydrodynamic interactions and the entanglements cannot longer be ignored. This is because at low concentrations the dilution is assumed to be infinite, therefore the polymer molecule is alone and has no interactions with other molecules.

4.3 Question 3

The biggest difference between these two viscometers is the third vent pipe at the Ubbelohde viscosimeters. This has the effect, that for the Ubbelohde-viscometer the effect of the sample is almost irrelevant. It also helps to get a constant hydrostatic pressure. Therefore the use cases vary a lot between those two. The Ostwald-viscometer is used for individual measurements and the Ubbelohde-viscometer is used for dialation measurements.

4.4 Question 4

The Staudinger-Index of the it-PMMA is higher because of the regular order of groups the molecule is sterically hindered. The at-PMMA has a random order therefore there is less steric hindrance. Another factor is the solvent because the it-PMMA is already more swollen the effect of the solvent is even stronger.